

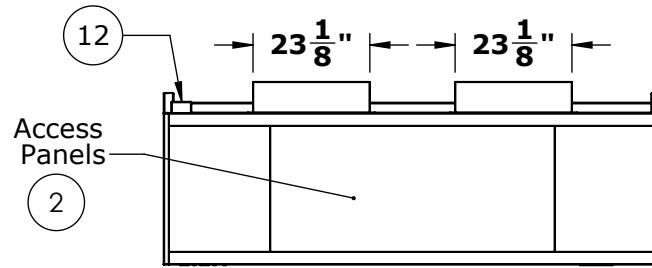
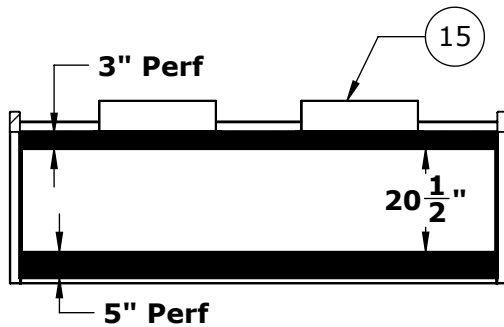
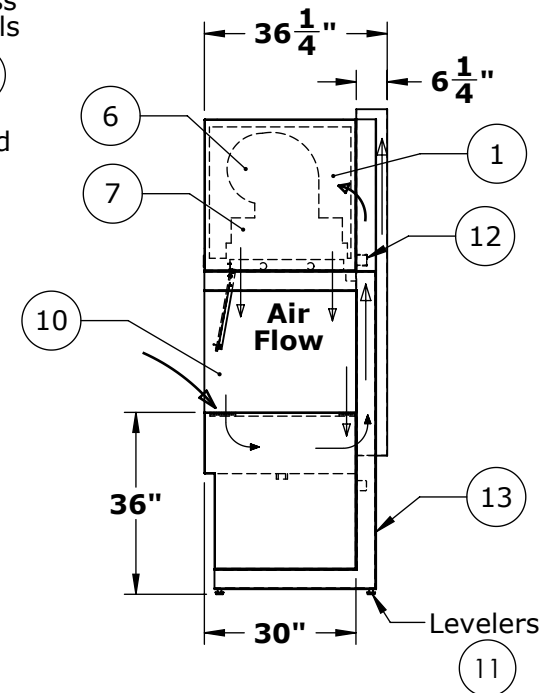



NOTES: UNLESS OTHERWISE SPECIFIED  
1. DIMENSIONS ARE IN INCHES (mm).

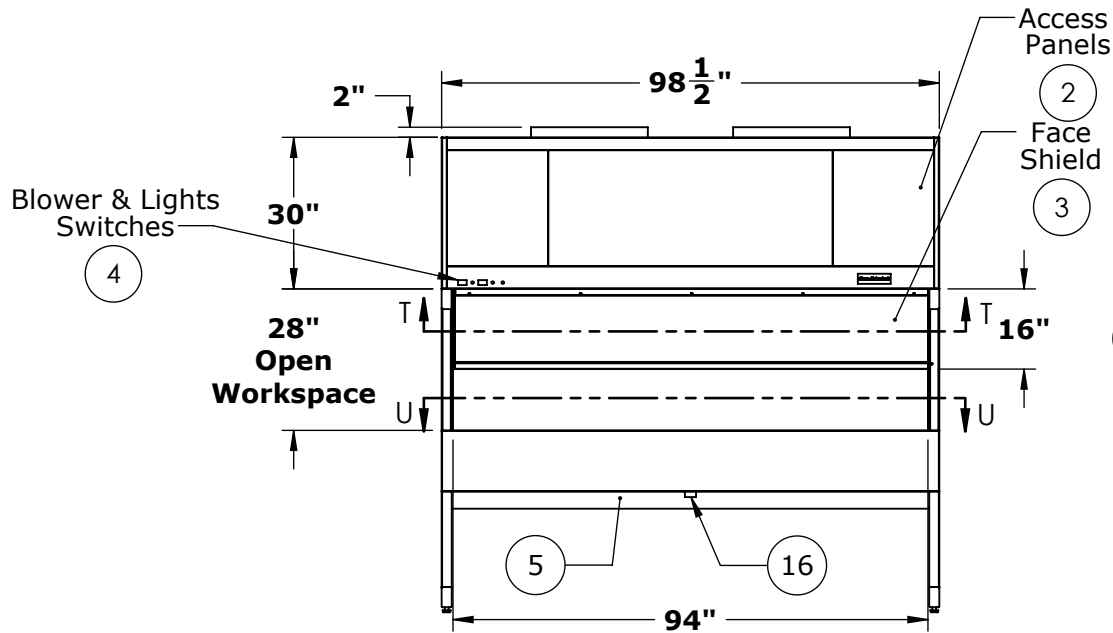

Top Cabinet ships separately  
from the lower frame

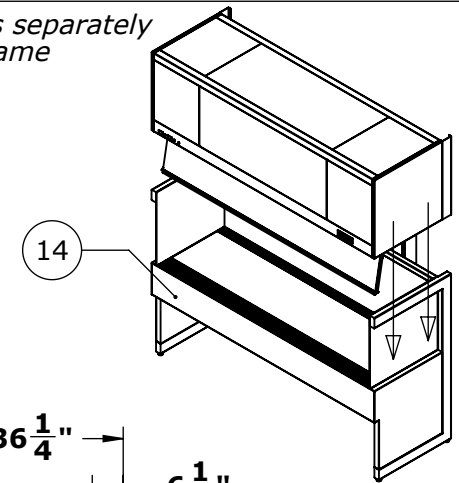

Work Surface  
SECTION U-U

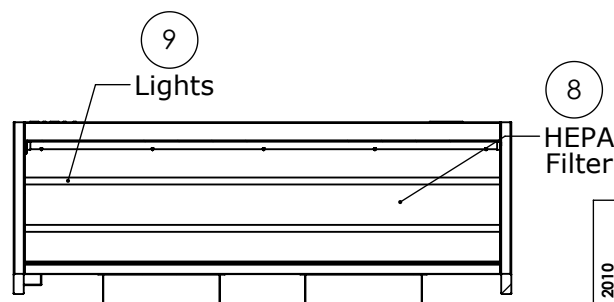

SECTION T-T



|                                        |                           |                                                                                                                                      |                      |          |                 |
|----------------------------------------|---------------------------|--------------------------------------------------------------------------------------------------------------------------------------|----------------------|----------|-----------------|
| © COPYRIGHT 2010<br>CLEAN AIR PRODUCTS | <b>Clean Air PRODUCTS</b> | 8605 WYOMING AVENUE NORTH • MINNEAPOLIS, MN 55445<br>800-423-9728 • ph: 763-425-9122 • fax: 763-425-2004<br>www.cleanairproducts.com |                      |          |                 |
|                                        | DESCRIPTION.              | <b>Vertical Flow Clean Bench</b><br><b>CAP412-8T30-36H-ER-SST</b>                                                                    |                      |          |                 |
|                                        | QUOTE NO.                 |                                                                                                                                      |                      |          |                 |
|                                        | <b>DO NOT SCALE</b>       | DATE DRAWN<br>09/07/2012                                                                                                             | QUOTE FORM<br>035505 | REV<br>1 | SHEET<br>1 of 2 |

NOTES: UNLESS OTHERWISE SPECIFIED

1. DIMENSIONS ARE IN INCHES (mm).

**Notes:**

1. Top Cabinet, 16ga CRS with White Powder Coat Finish
2. Blower and Filter Access Panels, Front and Rear  
(Field Personnel Can Relocate These In Order To Mount Prefilters In The Front)
3. Front Face Shield, Hinged 1/4" Polycarbonate, 12" Opening
4. Control Panel, Blower and Light Switch, Circuit Breakers
5. Back Panel, 18ga Stainless Steel
6. Double Blowers, 1/3HP, 120V, 60Hz
7. Supply Plenum, metal over top of HEPA, Flex duct to blower
8. HEPA, 99.99% on 0.3 micron particles, Aluminum anodized frame
9. Lights, Fluorescent Lamps, F96T12 (2x)
10. End Shields, Clear Polycarbonate
11. Leg Leveling Glides, Adjustable +/- 1/2"
12. Electrical Junction Box, (120V, 60Hz, 10.55 FLA)
13. Lower Support Frame, "C" Style, HRS with White Powder Coat
14. Stainless Steel Exhaust Plenum with 16ga Stainless Steel Work Surface
15. Exhaust Duct, Stainless Steel - Customer Supplied Exhaust must Provide a Minimum of 1200 CFM Total
16. Plenum Drain, 1-1/2" NPT Stainless Steel - Customer to Provide Plug if Needed

**Optional:**

17. Power Cord - 14 ft Long
18. 304 #4 Finish Stainless Steel Back Panel
19. 304 #4 Finish Stainless Steel Side Panel
20. 304 #4 Finish Stainless Steel "C"-Frame
21. 304 #4 Finish Stainless Steel Top Cabinet
22. Low Profile - 20" Tall Top Cabinet (Normally 30" Tall)
23. ULPA Filters - ISO Class 4, Fed Std 209E 10

**Shipping Note:**

1. Top Cabinet Shipped  
Separate from "C"-Frame

|                                        |                                                                                                                                                                |                      |          |                 |  |
|----------------------------------------|----------------------------------------------------------------------------------------------------------------------------------------------------------------|----------------------|----------|-----------------|--|
| © COPYRIGHT 2010<br>CLEAN AIR PRODUCTS | <b>Clean Air PRODUCTS</b> 8605 WYOMING AVENUE NORTH • MINNEAPOLIS, MN 55445<br>800-423-9728 • ph: 763-425-9122 • fax: 763-425-2004<br>www.cleanairproducts.com |                      |          |                 |  |
|                                        | DESCRIPTION. <b>Vertical Flow Clean Bench<br/>CAP412-8T30-36H-ER-SST</b>                                                                                       |                      |          |                 |  |
|                                        | QUOTE NO.                                                                                                                                                      |                      |          |                 |  |
|                                        | <b>DO NOT SCALE</b>                                                                                                                                            |                      |          |                 |  |
|                                        | DATE DRAWN<br>09/07/2012                                                                                                                                       | QUOTE FORM<br>035505 | REV<br>1 | SHEET<br>2 of 2 |  |